SPORT-ELEC®

I N S T I T U T

CELLU ROLL MASSAGER

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Réf : PM789

FR

GB

ES

NL

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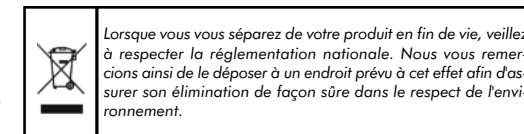
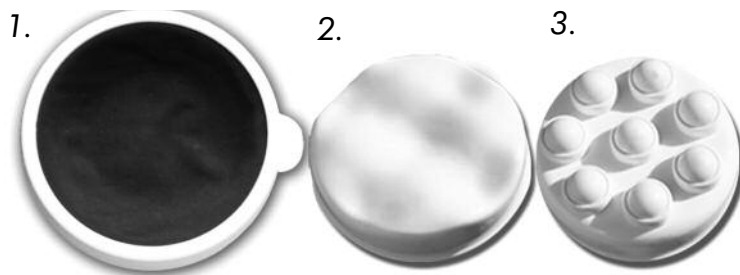
USER MANUAL

PRODUCT CARE AND MAINTENANCE

1. Place the CELLU ROLL MASSAGER in a safe cool and dry place when not in use.
2. Wipe with a damp soft cloth using no abrasive cleaner.
3. Never submerge the unit in any liquid.
4. Keep away from all solvents and harsh detergents.
5. DO NOT attempt to repair the CELLU ROLL MASSAGER. There are no user-serviceable parts.

SAVE THESE INSTRUCTIONS

Thank you for purchasing this CELLU ROLL MASSAGER. Before using this massager, please read the operation manual thoroughly.

CAUTION

- To reduce the risk of electric shock, do not remove cover or back. To reduce the risk of fire or electric shock, do not expose this product to rain or moisture.

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS BEFORE USING

All of the safety and operating instructions should be read, adhered to and followed before the product is operated.

TO REDUCE THE RISK OF ELECTRIC SHOCK, BURNS, FIRE OR INJURY :

1. Always unplug the appliance from the electrical outlet immediately after using and before cleaning.
2. Do not use while bathing or in a shower.
3. Do not place or store product where it can fall or be pulled into a tub or sink.
4. Do not place the product in the water or other liquid.
5. Do not reach for product that has fallen into water. Unplug it immediately.
6. Do not install roller where the hair may be easily twisted

WARNING

1. Do not leave product unattended when plugged in. Unplug it from outlet when not in use.
2. Do not leave product unattended when used by or near children.
3. Do not use attachments not recommended by manufacturer.
4. Never operate this product if it has a damaged cord or plug, if it does not work properly, if it has been dropped or damaged or dropped into the water.
5. Keep cord away from heated surfaces
6. Never block the air openings of the product or place it on a surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint and hair.
7. Never use while sleeping.
8. Never drop or insert objects into openings
9. Do not operate where aerosol products are being used or where oxygen is being administered.
10. Never drop or insert objects into openings.
11. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock or injury.
12. Do not carry this product by supply cord as handle.
13. Do not use on one concentrated area of the body. It is meant to be used in continuous motion, constantly changing the massage area.
14. To disconnect, turn all controls to the off position and then remove the plug from the outlet.
15. Do not perform any servicing on the product.
16. This product is designed for household use only.

CAUTION

1. Consult your doctor before using this product.
2. Do not use in the presence of unexplained pain.
3. Never leave the product operating or unattended, especially if children are present.
5. The appliance is designed as a non-professional massager to soothe tired and aching muscles. It is not intended to be a substitute for medical attention. If your ailments or symptoms persist or worsen after use, contact your doctor immediately.

INSTRUCTION FOR USE


- 1 - Make sure the power is off before plugging the device.
- 2 - Choose the accessory:
 - Roll-on massage head for anti-cellulite massage
 - Wavy massage head to relax muscles tensions
 - Protect net for the neck
 You can also use the roll-on or wavy massage head under the protect net.
- 3 - Apply cream or gel massage of your choice on your skin
- 4 - To turn on or off, turn the power button until "Click"
- 5 - To increase the intensity, slowly turn the power counterclockwise and clockwise to decrease the intensity, turn the power button in the direction of clockwise.
- 6 - Enter the aircraft through the wrist and start your massage with circular motions or going up and down the body.
- 7 - After use, remove the attachment and store your device in a dry and cool place.

Note: Unlock the safety net for the clean water if necessary. Dry before the screw on the device.

SPECIFICATION

Power: 220-240V
 Frequency: 50Hz
 Power: ≤ 15W
 Weight: 0.780 Kg

CE 2008 104 EC CEM

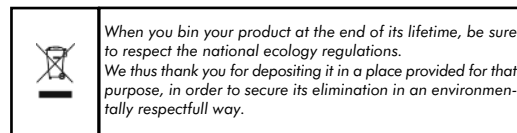ROHS 2002 95 EC This product is manufactured without materials which could have a harmful effect on ecology.

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SPECIFICATION

Power: 220-240V
Frequency: 50Hz
Power: ≤ 15W
Weight: 0.780Kg

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WARRANTY

This appliance is guaranteed for 12 months from the date of purchase. During the period of the guarantee, the distributor will take responsibility for defects occur, the appliance will be repaired or exchanged as the distributor deems necessary. Damage arising from improper use or normal wear and defects, which do not materially affect the performance of the appliance, is not covered under this guarantee. The guarantee is not valid if unauthorized repairs have been attempted or if replacement parts used are not approved by the supplier. The guarantee is valid only if it is stamped and signed by the retailer at the time of the purchase. For all claims made during the period of the guarantee, send the appliance and the warranty card to your retailer. All other claims are excluded from the guarantee unless specifically prohibited by law.

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